

AUTOSCALER®

THE SCIENCE OF DENTISTRY

 **South East Instruments**
A Pac-Dent Company

Trusted Quality & Performance

- Best-in-class warranty with 6 years for circuit boards, and 2 years on all other components
- Made for special periodontal and any supra-subgingival or routine scaling
- Linear power control allows exact power level needed for each procedure
- Automatic tuning circuit provides even power under all load conditions
- Compatible with most industry-leading magnetostrictive tips
- Tabletop models include hardwired handpiece assembly
- 25K and 30K models available in tabletop or built-in retrofit models

45
YEARS
of
ULTRASONIC EXPERIENCE



Order Information

Item#	Description
AW25C	Autoscaler® Ultrasonic Scaler, Tabletop Model, 25K
AW30C	Autoscaler® Ultrasonic Scaler, Tabletop Model, 30K
AW25BI	Autoscaler® Ultrasonic Scaler, Build-In Model, 25K
AW30BI	Autoscaler® Ultrasonic Scaler, Build-In Model, 30K

AUTOSCALER®

NEW!

Ultrasonic Inserts

 **South East Instruments**
A Pac-Dent Company



Universal (TFI-10)


Ultrasonic Insert Tip



Tip	Straight shank with medium length and tapered tip.
Deposit Types	• Light to heavy deposits • Tenacious calculus • Orthodontic cement
Applicable Areas	• Supragingival • Subgingival if tissue allows
Power Setting	Low to High

Thin

Ultrasonic Insert Tip

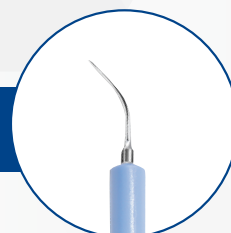

Tip	Straight shank with long length and a thin tapered tip.
Deposit Types	Light to moderate deposits
Applicable Areas	Subgingivally for periodontal pockets, but also can be used supragingivally
Power Setting	Low to Medium

Order Information

Item#	Description
INS 508601	Autoscaler® Universal Ultrasonic Insert Tip (TFI-10), 25K
INS 508701	Autoscaler® Universal Ultrasonic Insert Tip (TFI-10), 30K
INS 508602	Autoscaler® Thin Ultrasonic Insert Tip, 25K
INS 508702	Autoscaler® Thin Ultrasonic Insert Tip, 30K